POCKET AIR®


GOCare

GOCare Portable Nebulizer

- English
- German

EN

DE

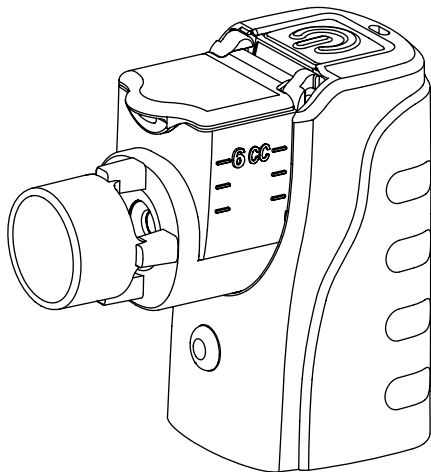Instruction Manual

MBPN002

CE
2797

Table of Contents

■ General Information	1
■ Intended Purpose	2
■ Important Safety Precautions	3
■ Explanation of Symbols	5
■ Package Contents	6
■ Optional Parts and Accessories	6
■ System Overview	7
■ How to Set Up the Medication Cup	8
■ Connecting to Power Supply	9
■ How to Inhale Properly	11
■ Cleaning and Disinfection	12
■ Troubleshooting	14
■ Specifications	15
■ Technical Data	16

General Information

EN

Thank you for choosing our GOCare Portable Nebulizer. GOCare is a high-quality and reliable handheld medical device designed to aerosolize medication for respiratory treatment.

MBPN002

The Aerosol is generated by active vibrating mesh technology which produces stable and consistent fine particles.

Not allowed to nebulize the medication that is not for respiratory illness prescribed by the doctor.

Benefits

- Increased lung function with FEV1 improvement
- Symptom severity score reduction
- Decreased healthcare utilization
- Treatment compliance increased

Please read this instruction carefully. The instructions for use are also available on the company website in multiple languages.

Please visit: microbasetech.com

Intended Purpose

Intended Use	The nebulizer is intended to turn the liquid (solution, suspension) into a mist for inhalation therapy. The device nebulizes the saline solution into aerosols for inhalation by the patient to soothe the discomfort, or prescription-only nebulization medicines. This device is not intended for connection to breathing equipment.
Indication for use	The indication for use of our device works to alleviate respiratory illness such as bronchitis or Chronic Obstructive Pulmonary Disease by inhalation therapy for long-term care use.
Contraindications	Patients who are unconscious, not breathing spontaneously or have pulmonary edema.
Intended user	This device is designed for all ages of patients suffering from respiratory illness such as bronchitis or Chronic Obstructive Pulmonary Disease for long-term care use. Except for patients who are unconscious, not breathing spontaneously or have pulmonary edema. The PATIENT is an intended OPERATOR. Children under 12 years old should be operated by an adult.
Environment	This device is intended for use in medical facility, such as hospital, clinic and doctor's office, and in general household. Temperature range: 10°C to 40°C Humidity: 30 to 85 % RH
Precautions for use	Please read the warnings and cautions in this manual carefully. And notice that any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/ or patient.
Validity period (Shelf Life)	Main unit 2 years Medication cup 1 year

Important Safety Precautions

READ ALL WARNINGS AND INSTRUCTIONS BEFORE USING THIS DEVICE.

GOCare Portable Nebulizer is a medical device. This section provides the information of warnings and cautions on using the product safely and correctly to prevent you and others from injury and device damage.

WARNINGS

EN

- Only use this device after consultation with a medical pharmacist or a doctor.
- Only use this device with medications prescribed by your doctor.
- The nebulizer is only intended for respiratory therapy, and any other application of this device is improper and dangerous.
- Keep the device out of unsupervised infant and children's reach.
- If any physical discomfort occurs, please stop using the nebulizer immediately.
- Do not share your nebulizer with others to avoid risk of infection. This product must only be used by a single patient.
- Do not keep product and accessories in a damp environment. Contamination and residual moisture encourage the growth of microorganism and increase the risk of infection.
- Close supervision should be provided by a parent or caregiver when this device is used by, on, or near children or disabled individuals to avoid other people inhale medicine and risk of injury and risk of strangulation with cable.
- For safety reasons, always disconnect the power adapter from the socket to avoid risk of electric shock under the following circumstances:
 - if a malfunction occurs during operation
 - before cleaning the device
 - immediately after use
- Be sure the device was properly cleaned before use to avoid possible contamination.
- Do not handle the power adapter with wet hands to avoid risk of electric shock.
- Do not swallow small parts of nebulizer to avoid risk of choking and injury.
- Cleaning and Disinfection
 - Clean and disinfect nebulizer according to instruction manual to avoid risk of infection.
 - Clean the medication cup after each use and disinfect the medication cup twice a week.
 - Allow the product's parts to dry after cleaning and disinfecting before reassembling the product and store in a clean location.
- Do not position the device so that it is difficult to disconnect from the supply mains (such as the mains plug or appliance coupler). Always ensure that the power source is easily accessible for quick disconnection.
- Use only accessories, detachable parts, and materials specified in this instruction manual.

Important Safety Precautions

CAUTIONS

- Do not attempt to clean the mesh with a cotton swab or any of reign objects, it may damage the mesh.
- Do not drop the device on the floor.
- Keep the device away from direct sunlight, excessive heat or cold to avoid damaging the batteries.
- Do not disassemble or attempt to repair the main unit while the ME EQUIPMENT is in use.
- Do not add more than 6 mL of medication to the medication cup.
- Do not start the nebulizer without filling in the medication.
- Do not poke the mesh with a cotton swab or a pin.
- Do not spill water or any liquid on the main unit, adapter or USB cable to avoid risk of electric shock.
- Do not pull the adapter out of the socket by the USB cable.
- Follow local laws and recycling plans regarding disposal or recycling of components, batteries and packaging.
- The device has only a start function, and neither the patient nor the operator can adjust any other functions.
- Allow the equipment to acclimate at 28°C ambient temperature for at least 1 hour after removal from minimum or maximum storage temperature before intended use.


Explanation of Symbols

Symbol / Meaning	
	Warning / Caution
	Class II equipment per IEC 60601-1
	Type BF equipment per IEC 60601-1
IP22	Protected against foreign objects equal to or greater than 12.5 mm in diameter and against drops of water falling at up to 15° from vertical
	Refer to instruction manual / booklet
	Do not re-use
	DC power source (direct current)
	Do not use if package is damaged
	Temperature limit
	Humidity limit
	Atmospheric pressure limitation
	Manufacturer
	Date of manufacture
	Serial number
	Batch code
	WEEE; waste electrical and electronic equipment. Do not treat this product as household waste.
	This device complies with the requirements of the Medical Devices Regulations (2017/745).
	Authorized representative in the European Community
	Medical device


Package Contents

Check Before Use

The following items are contained in the package.

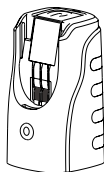1. Main Unit
MBPN00200

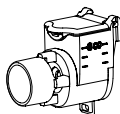2. Medication Cup
MBPN00201

Optional Parts and Accessories

Main Unit
MBPN00200

Medication Cup
MBPN00201

Mouthpiece
MBPN00207

(Single use)

Adult Mask
MBPN00203

(Single use)

Child Mask
MBPN00204

(Single use)

USB Power Adapter
MB0700403



USB Cable
MBPN00208

* If there are any missing parts, defects (e.g. sharp edges, rough surfaces), or malfunctions, please contact the store or distributor.

System Overview

EN

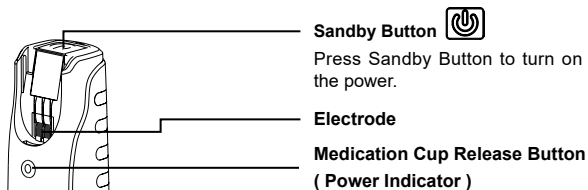
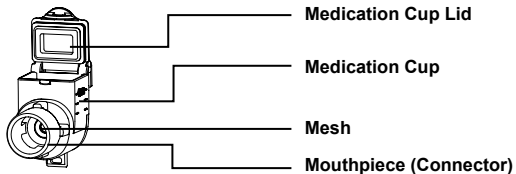
Main Unit

Press Release Button to remove the medication cup.

Status of Light Indicators		
 Light Indicators	Continuous Green Light	Nebulizer working normally
	Continuous Orange Light	Low batteries
	Flashing Orange Light	Malfunction
	No Light	No power, malfunction

▪ Bottom View

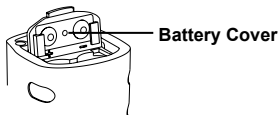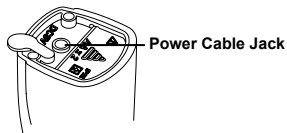How to Set Up the Medication Cup

Before using the product for the first time, read the instruction manual carefully.


Clean the medication cup and before using the product for the first time, after each use and after the product has not been used for an extended period of time and when parts are contaminated.

After adding the medication, avoid overturn the medication cup to prevent spillage.


1. Fill Medication

- Fill the medication in accordance with instruction from doctors.
The maximum capacity is 6mL.
- Make sure the medication cup lid had closed properly.
- Use a dry lint free cloth to wipe excessive medication or solution before install the medication cup onto the main unit.

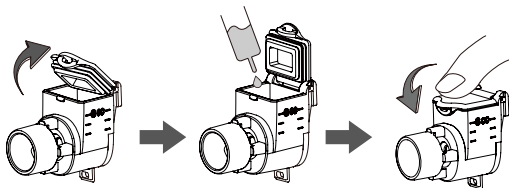2. Install the Medication Cup

- Align the buckle on the back of the medication cup with the trenches on the main unit and insert the medication cup. Check to ensure that the buckle is aligned and inserted correctly.

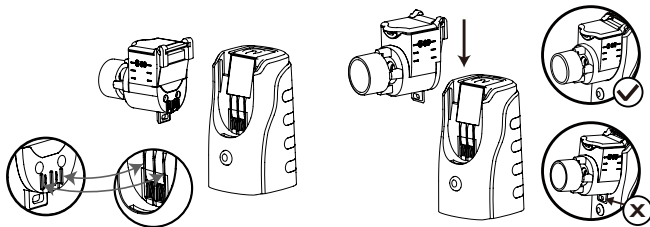
Connecting to Power Supply

Installing the Batteries

This device can be operated with two AA alkaline batteries.

EN

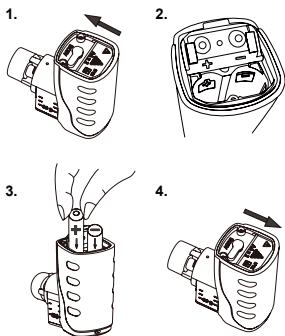- Press and push the battery cover in the direction of the arrow.
- Open the battery cover.
- Insert two AA alkaline batteries according to the battery polarity indicated inside the battery compartment.
- Push the battery cover in the direction of the arrow and push it until it catch the claws and clicks.

- 1. Do NOT insert a battery with the positive(+) and negative(-) poles reversed.**
- 2. Do NOT use old and new batteries together or mix and use different brand of batteries together.**

Battery life and battery replacement


1. When the power indicator turns orange (means low power), please replace both batteries with new one.
2. Battery life depends on the capacity, condition and storage environment condition of the batteries.

Connecting to Power Supply

Using the USB Cable with Adapter


1. Only use an adapter that is in compliance with the requirements of IEC 60601-1.
2. Please use the DC 5V/1A adapter, and the connection jack should be 9.5mm in length, with an outer diameter of 3.5mm, interior diameter of 1.35mm, and positive polarity.
3. The adapter does not charge batteries. As with all electronic devices, it is recommended to keep the product unplugged when not in use.
4. Please remove the batteries before using the adapter.
5. Keep the product unplugged when not in use.
6. Do not use speed charger.

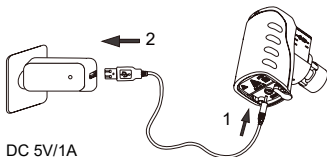
- Insert the jack connector of the USB cable into the cable jack on the main unit. Plug the USB cable into a standard adapter

How to Inhale Properly

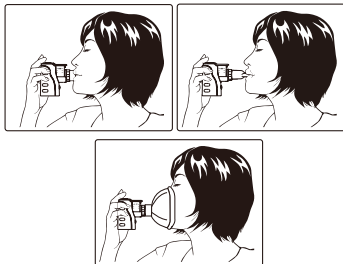1. Performing the inhalation

- Press Sandby Button on the top to start nebulization (light indicator turns green).
- Place the mouthpiece between your teeth, with your lips firmly sealed around the mouthpiece.
- * Accessories such as additional mouthpiece and mask could help to provide a more comfortable inhalation treatment.
- * If mask is used, gently press the mask against the face so that it fits snugly over the mouth and nose.
- Breathe in and out slowly through your mouth until aerosol formation stops.

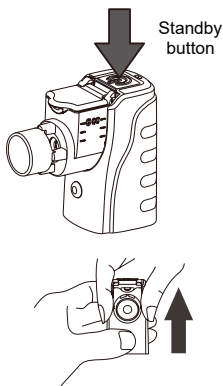2. End of Inhalation

- Press Sandby Button to turn it off.
- Press Release Button to remove the medication cup.
- Remove the batteries from the main unit.
- Discard the remaining medication and rinse the medication cup with distilled water and leave it air dry.

1. Tilt the device slightly until the solution is almost exhausted.
2. Do not re-use masks and mouthpiece to avoid risk of infection.

Cleaning and Disinfection

Device should be cleaned after every application to prevent the growth of pathogen, which increases the risk of infection.

After each use:

1. Remove the medication cup by pressing and holding the release button and pulling it upwards.
2. Remove the residual solutions and pour a small amount of distilled water to nebulize for 1-2 minutes to remove residuals completely.
3. Remove batteries or disconnect the adapter from the main unit.
4. Shake off excess water and allow parts to be fully air dried on a clean, dry towel.
5. Clean the main unit by wiping it off with a lint free cloth if necessary.
6. Make sure that all cleaned parts are completely dry before you store or use them next time.

1. Keep the battery compartment dry all the time.
2. Do NOT poke the mesh with finger, cotton swab or any object.
3. Do NOT clean parts in a dishwasher.
4. Do NOT use microwave to dry any parts.

Daily disinfections

Disinfection by alcohol


- It is important to disinfect the medication cup twice a week (disinfection once every 3 to 4 days).
- Rinse the medication cup with distilled water.
- Immerse the medication cup in 75% ethyl alcohol for one (1) minute.
- Rinse the medication cup with distilled water again, shake off excess water and allow parts to be fully air dry on a clean dry towel.
- Make sure that all parts are clean and dry before storage or use.

Cleaning and Disinfection


Alcohol is highly flammable. Do Not use alcohol near the fire or smoke to avoid risk of burn.

The following item can be disinfected by alcohol.

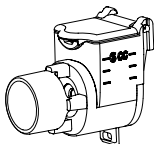Medication Cup

The following parts can be cleaned only with lint free cloth.

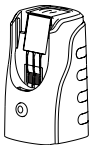• Main Unit


• USB Power Adapter

• USB Cable

Keep the main unit away from any liquid to avoid risk of electric shock.

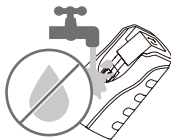Storage

Store the device and medication cup in a dry and clean environment.

1. Do NOT leave residual liquid in the medication cup, or re-use mask and mouthpiece to avoid risk of infection.
2. Remove the batteries if the device will not be in use for a prolonged period. Failure to do so may result in damage due to battery leakage.
3. Do NOT leave the device under direct sunlight, or in high humidity, extreme heat or cold environment.
4. Keep this device away from fire, high electromagnetic fields and out of reach of children.

Troubleshooting

In case of any of the below problems occur during use, please cross check with new Main unit and new Medication cup first, or refer to the table below:

Trouble	Possible Cause	Corrective Action
No/Low Nebulization	Poor power connection	Reconnect the adapter to the Main unit
		Confirm whether the negative electrode of the battery is sunken
		Cleaning the electrodes on Main unit with Rubbing alcohol for stain/rust remove
		Cleaning the electrodes on Battery cover with Rubbing alcohol for stain/rust remove
	Mesh clogging	Doing disinfection by alcohol, try again after Medication cup had completed dry
		Replace with new one
	The mesh is damaged	Replace with new one
	Aging of expired product	Replace with new one
No Lights, No Nebulization	The battery is not installed correctly	Follow (+) and (-) mark on the battery cover to re-insert the batteries in correct direction
	Batteries are depleted	Replace with new two AA alkaline batteries or connect the adapter instead
	Poor power connection	Reconnect the adapter to the Main unit
	Controller inside or battery compartment gets rusty	Cross confirm with new Main unit and new Medication cup
Green light on without nebulization	The mesh holes are clogged	Refer to Cleaning and Disinfection procedure to clean the medication cup
	The mesh is damaged	Replace the medication cup
Orange light on	Low batteries	Replace new two AA alkaline batteries or use the adapter instead
	Main unit malfunction	Replace with new Main unit
Orange light flashing	Poor power connection	Use alcohol pad to clean the electrodes and reinstall the medication cup
	Mesh deformed	Replace the medication cup
	Aging of expired product	Replace the medication cup
	The medication cup is not installed properly	Refer to How to Set up the medication cup procedure to reinstall the medication cup
	Turn on the power without Medication cup	Install the Medication cup before turn on the power

1. If the device does not nebulize normally after taking the above-mentioned procedure, please contact a healthcare professional, the nearest dealer or store where you purchased the device.
2. Never operate this device if any of the parts are not working properly or have been damaged to avoid risk of electric shock.
3. It is recommended to replace the medication cup after 1 year of use, and the device after 2 years of use.

Specifications

Product	GOCare Portable Nebulizer
Model	MBPN002
Method of operation	Active Vibrating Mesh Technology
Power supply	1.5 V AA Alkaline Battery x 2 5 V DC, 1 A USB Power Adapter*
Power consumption	< 1.5W (NOTE)
Vibrating frequency	117 kHz +/- 15%
Nebulization rate	≥ 0.25 mL/min
Particle size	MMAD < 5µm (NOTE)
Capacity of medication cup	Approx. 6 mL
Dimension	Approx. L78 X W41 X H73 mm
Weight	Approx. 74g (Excluding batteries)
Operating temperature and humidity	10 ~ 40°C, 30 ~ 85 % RH, 800 ~ 1060 hPa
Storage and delivery condition	-20 ~ 70°C, 20 ~ 75 % RH, 800 ~ 1060 hPa

* The optional adapter should be in compliance with the requirements of IEC 60601-1.

NOTE:

Tested with 0.9% normal saline, and under normal temperature of 23°C.

Technical Data

According to the Standard ISO 27427:2023 Anaesthetic and respiratory equipment — Nebulizing systems and components, Annex C “Test methods for aerosol output and aerosol output rate”, and Annex D “Test methods for particle sizing”.

Model	MBPN00201
Aerosol Output Rate (mL / min)	0.181
Aerosol Output (mL)	0.887
Residual Volume (mL)	0.033
Percentage of fill volume emitted per minute (%/min)	8.92
MMAD (Mass Median Aerodynamic Diameter) (µm)	4.637
Geometric Standard Deviation (GSD)	1.841
% > 5 µm (%)	46.64
% 2 µm to 5 µm (%)	46.87
% < 2 µm (%)	6.49
Noise	< 50 dB (approx. 1 meter)

- Test with Albuterol 0.1% (M/V) concentration in 0.9 % sodium chloride solution
- Performance may vary depending on the type of agent (in particular, a suspension and/ or high-viscosity solution), nebulizer used, and patient population. Please contact the manufacturer or drug supplier if you need additional information.

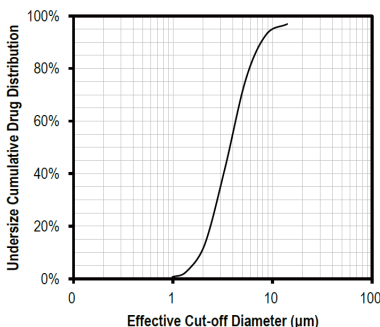
Figure1.

The performance of “Representative particle size distribution” for Albuterol as per ISO 27427:2023 Annex C and Annex D under adult breathing pattern.

Technical Data

In order to regulate the requirements for EMC with the aim to prevent unsafe product situations, the EN 60601-1-2 standard has been implemented. The nebulizer MBPN002 conforms to the EN 60601-1-2 standard for both immunity and emissions.

Nevertheless, do not use the nebulizer close to the strong electrical or electromagnetic fields. This may result in incorrect operation and create a potentially unsafe situation.

Guidance and manufacturer's declaration of MBPN002.

Manufacturer's declaration-electromagnetic emissions		
The <u>MBPN002</u> is intended for use in the electromagnetic environment (for home and professional healthcare) specified below.		
The customer or the user of the <u>MBPN002</u> should assure that it is used in such an environment.		
Emission test	Compliance	Electromagnetic environment-guidance (for home and professional healthcare environment)
RF emissions CISPR 11	Group 1	The <u>MBPN002</u> uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations /flicker emissions IEC 61000-3-3	Compliance	The <u>MBPN002</u> is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.

Technical Data

Manufacturer's declaration-electromagnetic immunity

The **MBPN002** is intended for use in the electromagnetic environment (for home and professional healthcare) specified below.

The customer or the user of the **MBPN002** should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment-guidance (for home and professional healthcare environment)
Electrostatic discharge(ESD) IEC 61000-4-2	Contact: ± 8 kV Air: ±2 kV, ±4 kV, ±8 kV, ±15 kV	Contact: ± 8 kV Air: ±2 kV, ±4 kV, ±8 kV, ±15 kV	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV for input/output lines	± 2kV for power supply lines Not applicable	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	± 0.5 kV, ±1 kV line (s) to line (s) ± 0.5 kV, ±1 kV, ± 2 kV line(s) to earth	± 0.5kV, ±1kV line (s) to line (s) Not applicable	Mains power quality should be that of a typical home and professional healthcare environment.
Voltage Dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	Voltage dips: 0 % UT; 0.5 cycle 0 % UT; 1 cycle 70 % UT; 25/30 cycles Voltage interruptions: 0 % UT; 250/300 cycle	Voltage dips: 0 % UT; 0.5 cycle 0 % UT; 1 cycle 70 % UT; 25 cycles Voltage interruptions: 0 % UT; 250 cycle	Mains power quality should be that of a typical home and professional healthcare environment. If the user of the MBPN002 requires continued operation during power mains interruptions, it is recommended that the MBPN002 be powered from an uninterruptible power supply or a battery.
Power frequency(50, 60 Hz) magnetic field IEC 61000-4-8	30 A/m 50 Hz or 60 Hz	30 A/m 50 Hz and 60 Hz	The MBPN002 power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

NOTE:

UT is the a.c. mains voltage prior to application of the test level.

Technical Data

EN

Manufacturer's declaration-electromagnetic immunity

The MBPN002 is intended for use in the electromagnetic environment (for home and professional healthcare) specified below.

The customer or the user of the MBPN002 should assure that is used in such and environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment-guidance (for home and professional healthcare environment)
Conducted RF IEC 61000-4-6	3 Vrms: 0.15 MHz – 80 MHz 6 Vrms: in ISM and amateur radio bands between 0.15 MHz and 80 MHz 80 % AM at 1 kHz	0.15 MHz – 80 MHz 6 Vrms: in ISM and amateur radio bands between 0.15 MHz and 80 MHz 80 % AM at 1 kHz	Portable and mobile RF communications equipment should be used no closer to any part of the MBPN002 including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance: $d = 1.2 \sqrt{P}$ $d = 1.2 \sqrt{P}$ 80 MHz to 800 MHz $d = 2.3 \sqrt{P}$ 800 MHz to 2.7 GHz Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m). Interference may occur in the vicinity of equipment marked with the following symbol: 
Radiated RF IEC 61000-4-3	3 V/m, 10 V/m 80 MHz – 2.7 GHz 80 % AM at 1 kHz	3 V/m, 10 V/m 80 MHz – 2.7 GHz 80 % AM at 1 kHz	
NOTE1 At 80 MHz and 800 MHz, the higher frequency range applies. NOTE2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

Technical Data

Recommended separation distance between portable and mobile RF communications equipment and the MBPN002

The MBPN002 is intended for use in an electromagnetic environment (for home and professional healthcare) in which radiated RF disturbances are controlled. The customer or the user of the MBPN002 can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the MBPN002 as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m		
	150 KHz to 80 MHz $d = 1.2 \sqrt{P}$	80 MHz to 800 MHz $d = 1.2 \sqrt{P}$	800 MHz to 2.7 GHz $d = 2.3 \sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where p is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

Technical Data

EN

Manufacturer's declaration-electromagnetic immunity

Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communications equipment

The MBPN002 is intended for use in the electromagnetic environment (for home and professional healthcare) specified below.

The customer or the user of the MBPN002 should assure that it is used in such an environment.

Test frequency (MHz)	Band ^a (MHz)	Service ^a	Modulation ^b	Maximum power (W)	Distance (m)	IMMUNITY TEST LEVEL (V/m)	Compliance LEVEL (V/m) (for home and professional healthcare)
385	380 – 390	TETRA 400	Pulse modulation ^b 18 Hz	1.8	0.3	27	27
450	430 – 470	GMRS 460, FRS 460	FM ^c ±5 kHz deviation 1 kHz sine	2	0.3	28	28
710	704 – 787	LTE Band 13, 17	Pulse modulation ^b 217 Hz	0.2	0.3	9	9
745							
780							
810	800 – 960	GSM 800/900, TETRA 800, IDEN 820, CDMA 850, LTE Band 5	Pulse modulation ^b 18 Hz	2	0.3	28	28
870							
930							
1,720	1,700 – 1,990	GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1, 3, 4, 25; UMTS	Pulse modulation ^b 217 Hz	2	0.3	28	28
1,845							
1,970							
2,450	2,400 – 2,570	Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, LTE Band 7	Pulse modulation ^b 217 Hz	2	0.3	28	28
5,240	5,100 – 5,800	WLAN 802.11 a/n	Pulse modulation ^b 217 Hz	0.2	0.3	9	9
5,500							
5,785							

Technical Data

NOTE:

If necessary to achieve the IMMUNITY TEST LEVEL, the distance between the transmitting antenna and the ME EQUIPMENT or ME SYSTEM may be reduced to 1 m. The 1 m test distance is permitted by IEC 61000-4-3.

- a. For some services, only the uplink frequencies are included.
- b. The carrier shall be modulated using a 50 % duty cycle square wave signal.
- c. As an alternative to FM modulation, 50 % pulse modulation at 18 Hz may be used because while it does not represent actual modulation, it would be worst case.

Technical Data

EN

Manufacturer's declaration-electromagnetic immunity Test specifications for ENCLOSURE PORT IMMUNITY to proximity magnetic fields

The MBPN002 is intended for use in the electromagnetic environment (for home and professional healthcare) specified below.

The customer or the user of the MBPN002 should assure that it is used in such an environment.

Frequencies	Test Level [A/m]	Modulation	Dwell time [s]	Compliance LEVEL [A/m] (for home and professional healthcare)
30 kHz ^a	8	CW	3	8
134.2 kHz	65	Pulse modulation ^b 2.1 kHz	3	65 ^c
13.56 MHz	7.5	Pulse modulation ^b 50 kHz	3	7.5 ^c

NOTE:

a. This test is applicable only to ME EQUIPMENT and ME SYSTEMS intended for use in the HOME AND PROFESSIONAL HEALTHCARE ENVIRONMENT.

b. The carrier shall be modulated using a 50 % duty cycle square wave signal. c.r.m.s., before modulation is applied.

WEEE (Directive on Wasted Electrical and Electronic Equipment)

This marking shown on the product indicates that it should not be disposed of, with other household wastes at the end of its working life. To prevent possible harm to the environment or human health from uncontrolled waste disposal, please separate this from other types of wastes and recycle it responsibly to promote the sustainable reuse of material resources.

Household users should contact either the retailer where they purchased this product, or their local government office, for details of where and how they can take this equipment for safe recycling.

MicroBase Technology Corp.

No. 756, 758 Jiadong Rd., Bade Dist., Taoyuan City, 33464, Taiwan.

Warranty card

Product Name	GOCare Portable Nebulizer
Model Number	
Serial Number	
Purchase Date	
Dealer	

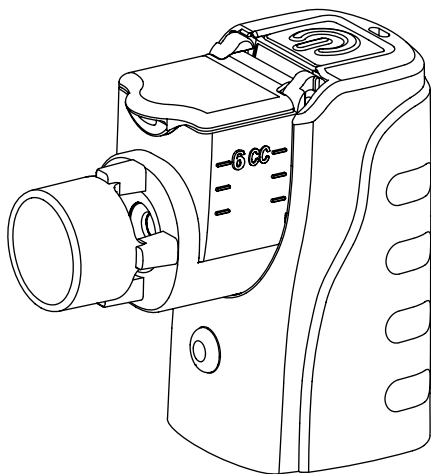
1. Please get a shop stamp from your dealer, put the purchase date on this card, and keep this card with you. Presentation of this card is required when a warranty service is requested. (No new card will be issued to you if this card is lost.)
2. MicroBase guarantees that the main unit of POCKET AIR, if used according to the instructions, will be free from defects in material and workmanship caused by the manufacturing process for the warranty period of 2 years. Any other replacement parts or accessories are not covered hereby. You will void this warranty if you disassemble the main unit of this GOCare Portable Nebulizer.
3. The warranty period of the medication cup will be 12 months, starting from the purchase date. The warranty does not apply to (what's NOT covered) :
 - I . Mesh clogging.
 - II . Any failure to follow our written instructions for the product.
 - III . Damage or defect due to willful neglect or negligence by anyone other than us.
4. During the warranty period, Microbase will repair or replace the device if a defect is discovered. If it is replaced, the replacement device may either be the same model or a model that is at least comparably equipped. Replacement or repair of the device shall not serve as the basis for a new warranty. All replaced old devices or parts shall become the property of Microbase. Further claims are excluded. This disclaimer of warranty shall be ineffective in the event of injury to life, limb or health, in cases of wrongdoing and gross negligence, product liability and if substantive obligations under the warranty agreement are violated. Any expenses of non-maintenance work (including but not limited to, the back and forth shipping cost or expenses arising from your specific request) will be borne by you.
5. The warranty shall be cancelled if:
 - the device has been operated or used improperly with respect to the descriptions in the instructions for use
 - damage is present that is attributable to the effects of water, fire, lightning, etc.
 - the damage was caused by transporting the device incorrectly or a falling impact
 - the device has been misused or not cared for correctly
 - repairs, adaptations or modifications have been made to the device by persons not authorized by Microbase
6. For any service requested without this warranty card, beyond the warranty period, or not covered by this warranty, a service fee will be charged.
7. In the event that replacement parts or warranty service is needed, please contact the nearest dealer or the store where you purchased this device.

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